

AUDIOMETERS

Cello

Diagnostic Audiometer




Orchestral perfection
A full line of Audiology Equipment

- Pure tone and speech audiometry in AC, BC, FF
- Video VRA option available
- High-frequency audiometry (optional)
- Wall-mounting option and dedicated transducers support available

DESIGNED, ENGINEERED AND MANUFACTURED IN ITALY



COMPUTER CONTROL

Connect Cello to the USB port of your Windows® computer and you are ready to go: no driver and no additional cable required. You can control Cello through the Maestro's amazingly user-friendly interface, available also for Noah.

TECHNICAL DATA

AUDIOMETRY

- Intensity: AC: -10 to 120 dB HL
BC: -10 to 80 dB HL
- Freq. range: 125-8000 Hz
8-20 kHz (opt.)
- Input: Tone, warble tone, EXT1, EXT2, USB audio, mic
- Masking: NBN, WN, SN
- Output: AC L+R, insert L+R, BC, FF L+R,
insert masking, HF L+R
- Talk over: External or computer microphone
- Monitor: External earphone/speaker or computer speaker
- Test battery: PTA (HL and UCL), speech audiometry (WRS, SRT, MCL, and UCL), autothreshold, SISI, Stenger, Master Hearing Aid, QuickSIN™ (opt.), High Frequency (opt.), video-VRA (opt.)

DEVICE SPECIFICATIONS

- Conn. to PC: USB
- Compatible sw: Inventis Maestro (stand-alone or Noah module)
- Controlling App: Maestro
- Dimensions: 30,8 x 21,7 x 6,4 cm
12.1 x 8.5 x 2.5 in. (DxWxH)
- Weight: 1.60 kg / 3.53 lbs


Cello is a PC-controlled diagnostic audiometer, capable of performing fast and accurate pure-tone and speech audiometry exams, through air, bone conduction, and free-field. It features several additional tests, including QuickSIN™, HF audiometry, and video-VRA. With its revolutionary design, complete feature set and top reliability, Cello is a little masterpiece, ideal for anyone unwilling to compromise on appearance and user experience.

AVAILABLE TESTS

Beside the standard pure-tone and speech audiometry, Cello includes Hughson-Westlake automated audiometry, high-frequency audiometry, QuickSIN®, SISI, Stenger and includes a 2 channel Master Hearing Aid.

VIDEO-VRA

With Cello, in few minutes you can also build a video-VRA system with unlimited possibilities for customizing the images and videos used as reinforces. Everything is controlled through the Maestro module, simply connect up to three displays to your computer to provide the reinforces.

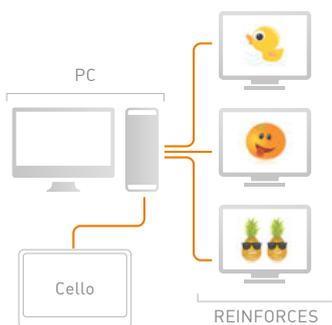

Find our distributors at:
www.inventis.it/world/dealers-inventis

Discover the complete orchestra at: www.inventis.it

